





VELA Tango

VERSION: 700 & 700E

Manual no. 105675

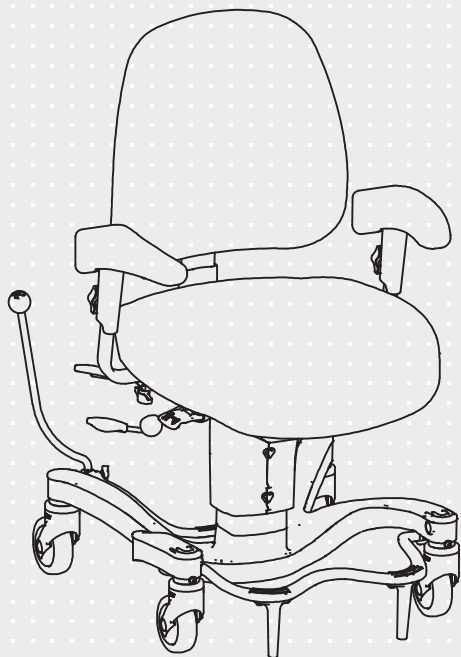

VELA Tango 700

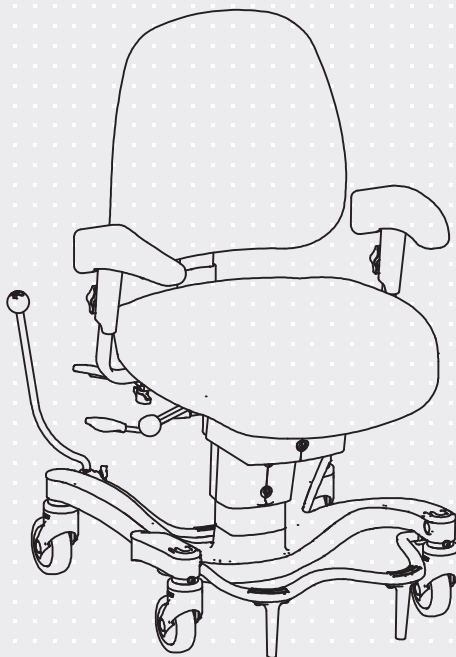

VELA Tango 700E



ISO 9001
ISO 14001

CHAIR FUNCTIONS



Backrest height



Backrest angle



Seat depth



Brake



Seat angle



Tilt



Armrest width




Armrest height



Seat rotation
(handle on left side)



Electric height adjustment

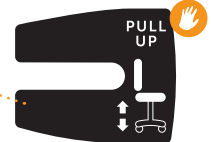

Manual height adjustment



GETTING STARTED

- 1 Brake the chair and sit back in the seat.
- 2 There should be about 3 cm between the seat cushion and the back of your knees.
- 3 Adjust the seat angle and height so that you sit comfortably with your feet flat on the floor.
- 4 Adjust the backrest height and angle to adjust the seat depth and support your lower back and back.
- 5 Adjust the armrest width and height so that your forearms rest on the armrest pad.

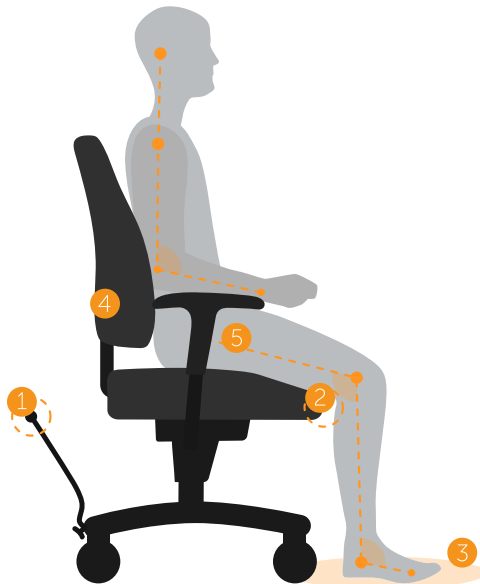


TABLE OF CONTENTS

1.0. INTRODUCTION	5
1.1. WARNINGS	6
1.2. TROUBLESHOOTING	7
1.3. VELA TANGO WITH ELECTRIC HEIGHT ADJUSTMENT	8
1.3.1. Power input	8
1.3.2. Emission	8
1.3.3. Surge protection	8
1.4. LABELS USED ON THE CHAIR	9
1.5. SAFETY	9
1.5.1. Intended use	9
1.5.2. Transport	10
1.6. WARRANTY	10
1.7. UNPACKING THE CHAIR	10
1.8. MAINTENANCE	10
1.8.1. Wheels	11
1.8.2. Frame	11
1.8.3. Padding	11
1.8.4. Recycling guide	11
1.9. RECYCLING AND DISPOSAL	11
2.0. SETTINGS	12
2.0.1. Attaching the backrest	12
2.0.2. Height adjustment (Electric)	12
2.0.3. Changing place for height adjustment*	12
2.0.4. Height adjustment (Manual)	13
2.0.5. Setting angle of seat	13
2.0.6. Adjustment of seat tilt	13
2.0.7. Setting angle of backrest	13
2.0.8. Adjustment of seat depth	14
2.0.9. Setting the back height	14
2.0.10. Back bolt*	14
2.0.11. Adjusting lumbar support*	14
2.0.12. Adjusting the armrest height	15
2.0.13. Adjusting the armrest width	15
2.0.14. Fold-down armrest*	15
2.0.15. Rotation	15
2.0.16. Brake	16
2.0.17. Foot pedal for braking	16
2.0.18. Setting the angle of the brake	16
2.0.19. Changing the seat	16
2.0.20. Charging (el)	17
2.0.21. Soft touch contact (El)*	17
2.0.22. Hand box (El)*	17
2.0.23. Electric brake (el)*	17
2.0.24. Folding footrest*	18
2.0.25. Side-mounted swivel footrest*	18
2.0.26. Strolling bracket*	18
2.0.27. 4-wheel brake*	18
2.0.28. Accessories	19

1.0. INTRODUCTION

DEAR CUSTOMER

Congratulations on your new VELA chair! We are sure this chair will give you daily pleasure. This user manual provides useful information on adjusting, operating and maintaining your chair.

Important

Please read this user manual carefully and keep it safe for when you need to refer to it. This manual is also available on VELA's website: www.vela.eu, where you can also find other relevant information.

If you have questions, experience errors or adverse events, please contact your VELA supplier.

We reserve the right to make changes.

VELA

Gøteborgevej 8-12
9200 Aalborg SV
Denmark
www.vela.eu

1.1. WARNINGS



Use	In case of an allergic reaction, stop using the chair and consult a doctor
Use	Only use original spare parts from VELA
Use	Only use the supplied charger for charging the chair
Use	Do not make any changes to the accessories that come with the chair
Use	Only use the chair for the intended purpose
Use	Only adjust the seat angle and backrest angle while the user is sitting in the chair and the user's feet are flat on the floor
Use	Always have the chair brake on when you sit down or get up
Use	The brake rod may only be used to operate the brake. If putting too much load on the brake rod, the brake rod/holder may brake and thereby lose its function
Use	On smooth surfaces such as highly polished floors or marble tiles, the chair may slide. Here, it is recommended to use a 4-wheel brake
Use	Never use the chair on surfaces with level differences as the chair may become unstable
Use	Maintenance other than cleaning should always be performed by a professional
Use	Never service the chair while it is in use
Use	Always follow the recommended maintenance and cleaning instructions
Use	The chair's functions are intended for the user. Only the user or trained professionals should adjust the chair while the user is sitting in it
Use	Never use the electric height adjustment for more than two minutes at a time
Use	Risk of shearing or squeezing fingers in moving parts. Therefore, never use the chair if someone is touching underneath the seat
Choking hazard	Choking hazard on small parts. Never allow children to play with the chair, as small parts are removable
Choking hazard	If the chair is equipped with electric height adjustment, cables will be present. Never let children play with these
Electric noise	Using cables other than the standard ones may increase electrical emission or decrease electrical immunity
Electric noise	Avoid placing the chair close to appliances that emit electrical noise
Electric noise	The chair may only be used with accessories mentioned in this user manual, as electromagnetic interference may otherwise occur
Electric noise	If the chair moves unintentionally or does not work as expected, find another location for the chair

1.2. TROUBLESHOOTING

Mistake	Possible cause and solution
Electric height adjustment does not work	<p>The plug is not in the armrest – plug it in</p> <p>There is no power on the chair's battery. Charge the battery with the supplied charger</p> <p>The fuse is broken. See section 1.3.3</p> <p>Perform a reset of the actuator. Hold both buttons (up and down) simultaneously for 20 seconds. Keep holding the buttons until the actuator has fully retracted and then moved approximately 5 mm up again.</p>
Manual height adjustment does not work	<p>There is too much load or uneven load on the chair. Try to adjust the height without sitting in the chair. If this does not work, the gas cartridge is defective and should be replaced. Contact your VELA supplier</p>
The brake function does not work	<p>The brake is not activated. See section 2.0.16-2.0.18., 2.0.23. or 2.0.27.</p>
The armrest seems loose	<p>One of the hand screws might be loose. Tighten the hand screw</p>
The chair rolls badly	<p>Wheels are blocked. Clean the wheels for wire residue and the like.</p>
Brake rod or chassis is in the way when I move in the chair	<p>You have a chair with seat rotation. Use seat rotation to straighten the chair. See section 2.0.15</p>
Brake rod seems loose	<p>The hand screw holding the brake rod might be loose. Tighten the screw. See section 2.0.18.</p>

If the above does not solve your problem, you should contact your VELA supplier.

1.3. VELA TANGO WITH ELECTRIC HEIGHT ADJUSTMENT

1.3.1. POWER INPUT

Voltage: 24Vdc

Current: 0.5A

Power: 12VA

1.3.2. EMISSION

A VELA chair with electrical functions has been tested and complies with current requirements regarding electromagnetic radiation (EMC requirements). The chair is tested to the following specifications:

:: Emission class CISPR 11, Group 1, Class B

See more details at www.vela.eu

1.3.3. SURGE PROTECTION

A fuse is installed on the chair which secures the chair against power surges. The replaceable part of the fuse has the following specifications:

:: Flat fuse, type: ATOF blade fuse, fast-acting

:: Voltage: 32V




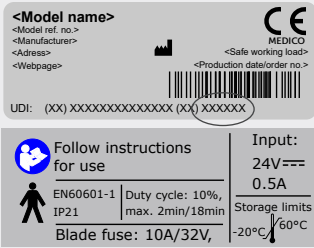






:: Current: 10A

:: Max. current: 1000A

Replacement of the fuse may only be performed by professionals, as the cause of the breach must be documented for safety reasons.

IN GENERAL

1.4. LABELS USED ON THE CHAIR

Warning	Risk of shear and squeeze: Keep hands and fingers away from moving parts						
Icon	Read "Instruction For Use" before use						
Icon	Part used: Type B cf. IEC60601-1 covering the following areas of the chair:						
Labels	60 x 22.5 mm The last six digits are the product ID number: This label only applies to chairs with electrical functions: The chair is IP21 approved, which corresponds to protection against water droplets as well as against touching dangerous parts with a finger.	 <p> <Model name> <Model ref. no.> <Manufacturer> <Address> <Webpage> </p> <p> CE MEDICO <Safe working load> <Production date/order no.> </p> <p> UDI: (xx) xxxxxxxxxxxxxxxxx (xx) xxxxxxxx </p> <table border="1"> <tr> <td> Follow instructions for use</td> <td rowspan="2"> Input: 24V= 0.5A </td> </tr> <tr> <td>  EN60601-1 IP21 Duty cycle: 10%, max. 2min/18min Blade fuse: 10A/32V, </td> </tr> <tr> <td colspan="2"> Storage limits -20°C / 60°C </td> </tr> </table>	 Follow instructions for use	Input: 24V= 0.5A	 EN60601-1 IP21 Duty cycle: 10%, max. 2min/18min Blade fuse: 10A/32V,	Storage limits -20°C / 60°C	
 Follow instructions for use	Input: 24V= 0.5A						
 EN60601-1 IP21 Duty cycle: 10%, max. 2min/18min Blade fuse: 10A/32V,							
Storage limits -20°C / 60°C							

Labels on the chair should be read at a distance of approx. 40 cm under normal daylight conditions.

For labels related to the operation of the chair, please refer to section 2.0.

1.5. SAFETY

This chair is CE marked and complies with applicable EU regulations. The chair is manufactured by: VELA, Gøteborgvej 8-12, 9200 Aalborg SV, Denmark.

The VELA Tango with electrical functions meets all safety requirements according to IEC 60601-1, including requirements for EMC according to IEC60601-1-2.

1.5.1. INTENDED USE

This device is intended to support a user in a sitting posture while performing indoor activities of daily living including work activities.

- :: The chair is not built to stand on
- :: The wheel brake must always be activated when sitting down or standing up

All variants of this chair are tested and approved for a maximum load of 160 kg. For optimal use, the chair's recommended user weight is 50-130 kg with gas cartridge and 50-145 kg with electric height.

Note: Persons using a harness in the chair must always be supervised when the harness is in use.

1.5.2. TRANSPORT

When transporting the chair by car, airplane or any other means of transport, always use the chair brake. The chair is not approved for use as a passenger seat in any means of transportation and may not be used as such.

The chair is designed to be transported under the following environmental conditions:

Temperature: - 20°C to + 60°C

Humidity: 15-95%

Allow the chair to stand for 24 hours before use if it has been stored below - 10°C.

1.6. WARRANTY

VELA grants a right of claim under applicable law in the country where the chair is acquired. The right of complaint applies only when using original spare parts and accessories, as well as customisations made by professionals.

Note: if an exclusive contractual agreement has been made directly with VELA, this applies!

VELA assumes no liability for damages to the product or the user caused by:

- :: Transportation
- :: Misuse
- :: Normal wear and tear
- :: Wrong use
- :: Use of spare parts and accessories that are not original VELA products
- :: Adjustments made by unauthorised service technicians

1.7. UNPACKING THE CHAIR

The chair comes in environmentally friendly and recyclable packaging.

In case of defects or transport damage to the chair, please contact your VELA supplier right away. In such cases, do not attempt to use the chair or to fix it yourself.

On delivery you have received:

- :: VELA chair
- :: Instructions for use
- :: Charger (only for chairs with electrical functions)

Note: see item 2.0.1 for correct attachment of backrest.

1.8. MAINTENANCE

We recommend an annual inspection of the chair to identify any defective or worn parts. All removable parts should be tightened regularly.

You can find various guides on www.vela.eu.

1.8.1. WHEELS

We recommend you regularly check the wheels and remove any residual threads, etc. Certain floor types accumulate dirt on the castors. Therefore, clean the castors as needed.

1.8.2. FRAME

Clean with lukewarm water and if necessary, add some pH-neutral detergent.

1.8.3. PADDING

Fabric: Clean by chemical cleaning, with furniture cleaner or with water using a hard-wrung cloth.

Synthetic fabric or leather: Clean with water or mild soap using a hard-wrung cloth.

Leather: Clean by vacuuming with a soft brush. Avoid using detergent.

Note: On www.vela.dk you can find cleaning instructions for all our standard upholstery types.

1.8.4. RECYCLING GUIDE

Your VELA chair can be reused by another person when you no longer need it as long as the chair is serviced by a technician beforehand.

The chair's life expectancy is five years for regular and proper use. We recommend that VELA, your VELA supplier or a trained professional perform an annual inspection with the replacement of worn parts.

1.9. RECYCLING AND DISPOSAL

This product contains recyclable materials. Therefore, the product must be disposed of in accordance with local regulations and not together with ordinary household waste. Proper disposal and recycling will help prevent negative effects on the environment and people.

Note: Electrical and electronic equipment and accompanying batteries contain materials, components and substances that can be harmful to human health and the environment if the waste is not handled properly.

Electrical and electronic equipment and batteries are marked with the icon shown below. This symbolises that electrical and electronic equipment and batteries must not be disposed of together with unsorted household waste but must be collected separately.



It is important that you return your used batteries to the established schemes. In this way, you help ensure that batteries are recycled in accordance with current legislation and do not cause unnecessary harm to the environment.



2.0. SETTINGS

The gas spring and the electric lifting mechanism should only be replaced by a VELA technician or a person trained in this field.




2.0.1. ATTACHING THE BACKREST

The backrest is lowered into the hole. Screw the handle back on and tighten securely.

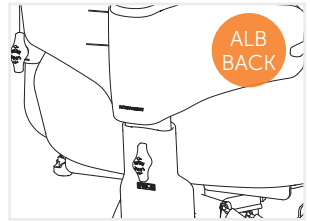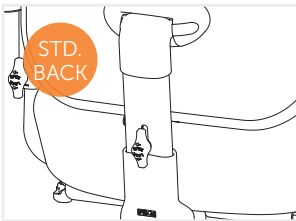

2.0.2. HEIGHT ADJUSTMENT (ELECTRIC)

Can be operated from the underside of either the right or left armrest. VELA Tango El is also available with height adjustment via hand box (if caregivers are to be able to adjust the height). **Note:** Height function should only be used continuously for one minute, then a nine-minute break should follow.

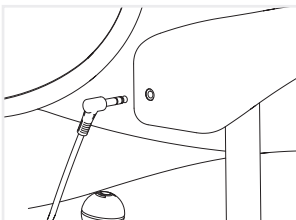

2.0.3. CHANGING PLACE FOR HEIGHT ADJUSTMENT*


Remove the plug from the back of the armrest and insert it into the opposite armrest or bottom of the hand box. **Note:** Height function may only be used continuously for one minute, then nine minutes of break.

Note: Setting the angle of backrest and seat should always be done with your feet flat on the floor.




2.0.4. HEIGHT ADJUSTMENT (MANUAL)

Adjusted by pulling up the front handle on the right-hand side.

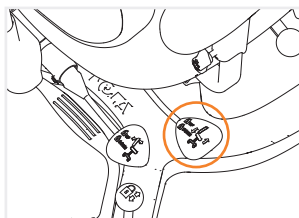


2.0.5. SETTING ANGLE OF SEAT

Adjusted by pulling up the middle handle on the right-hand side.

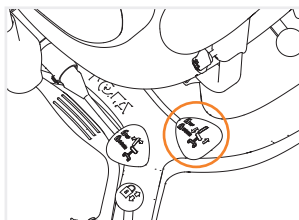

2.0.6. ADJUSTMENT OF SEAT TILT


The entire seating unit moves (seat and backrest together) by lifting the middle handle on the right side.

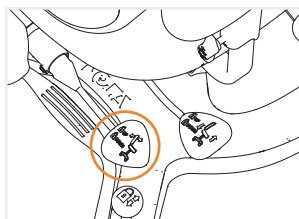

2.0.7. SETTING ANGLE OF BACKREST

Adjusted by pulling up the rear handle on the right-hand side. Only operate the handle when someone is sitting in the chair - as the backrest tilts forward.


* Depending on model

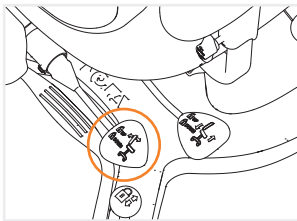

2.0.8. ADJUSTMENT OF SEAT DEPTH


Adjusted by pulling up the rear handle on the right side (backrest angle). Adjusted so that the seat depth suits the user of the chair.

Operate the handle only when someone is sitting in the chair, as the backrest tilts forward.

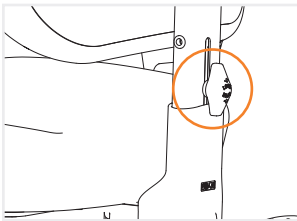


2.0.9. SETTING THE BACK HEIGHT

Release the handle located on the back. Adjust the back. Remember to tighten the handle afterwards.

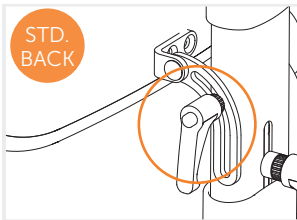


2.0.10. BACK BOLT*

The backrest can be set to a specific angle by adjusting and attaching the bolt mounted on the back of the backrest.

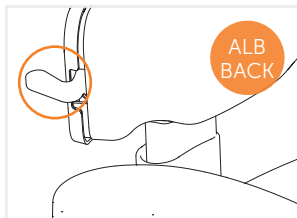


2.0.11. ADJUSTING LUMBAR SUPPORT*

The backrest can be adjusted to provide optimal lumbar support by pulling up the handle on the right-hand side of the backrest. The backrest adjusts to fit and support the lower back. Release the handle, and the backrest remains in the desired lumbar angle.

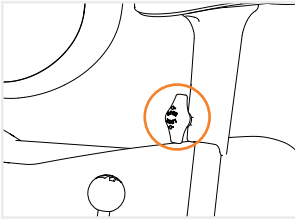


2.0.12. ADJUSTING THE ARMREST HEIGHT

Adjusted by loosening the screw on the armrest bar and adjusting the height of the armrest. Remember to tighten the screw when done.

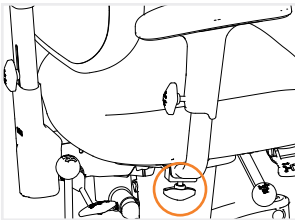

2.0.13. ADJUSTING THE ARMREST WIDTH


Adjusted by loosening the screw under the seat on the seat support and adjusting the armrest width by pulling out the armrest. Remember to tighten the screw when done.

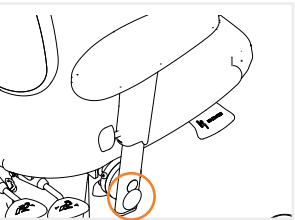

2.0.14. FOLD-DOWN ARMREST*

To fold down the armrest, press the button on the side of the armrest and push the armrest backwards and down.


Note: The armrest is only locked when it is in an upright position. The armrest can only be folded down if the width adjustment of the armrest allows it.

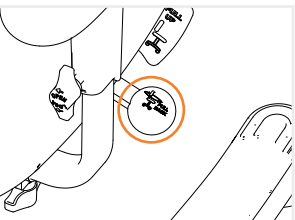

2.0.15. ROTATION

The seat is rotated by pushing back on the front handle on the left-hand side. Rotation can be turned and locked in the following positions: 90° to the left, 0° (neutral), and 90° to the right.


* Depending on model

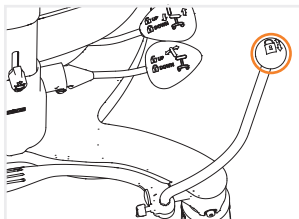

2.0.16. BRAKE


The brake is activated by moving the handle on the side of the chair forward. Unlock the chair by moving the handle back into the upright position.

Note: See Setting the angle of the brake, section 2.0.18. for adjusting the position of the brake and the distance to the seat.

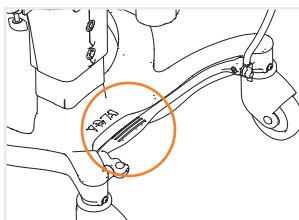

2.0.17. FOOT PEDAL FOR BRAKING


The brake can be activated with the foot pedal by pressing the foot against the brake shaft in the area shown.

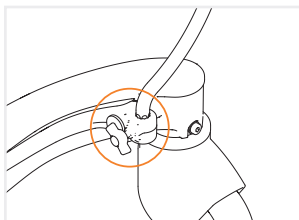

2.0.18. SETTING THE ANGLE OF THE BRAKE


Loosen the thumbscrew and remove the brake rod from the holder. Insert the brake rod in the desired side and angle. Remember to retighten the thumbscrew.

Note: The brake can be set at 30° intervals. The brake rod can be moved forwards/backwards and closer to/further from the seat as required. The brake can also be moved to the opposite side.

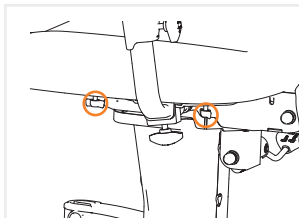

2.0.19. CHANGING THE SEAT

The seat can be replaced by removing the four thumbscrews that are located in each corner under the seat.

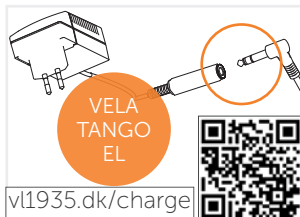

2.0.20. CHARGING (EL)

Charge every day. Remove the operation plug from the back of the armrest and place it in the charger as shown. The charger is plugged into an electrical outlet. When the chair is charging, the LED on the charger is yellow. When the chair is fully charged, the LED turns green.

Note: The chair must not be used during charging



2.0.21. SOFT TOUCH CONTACT (EL)*

Placed on the underside of the left or right-hand armrest as required.



2.0.22. HAND BOX (EL)*

The hand box is used via push buttons on the front. VELA hand boxes can have one or more functions. The icons on the hand box show the button functions.

Note: Insert plug in hand box to operate functions.




2.0.23. ELECTRIC BRAKE (EL)*

The electric brake is activated and deactivated using the hand box or switches. It is advisable to activate or deactivate the brake completely when using this feature. Can also be used together with the 4-wheel brake (see section 2.0.27.)


Note: The brake can only be used via switches under the armrest, the hand box, or the Soft touch contact.

* Depending on model

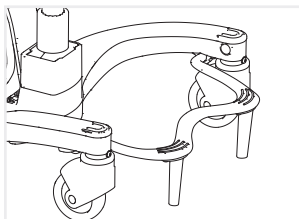


2.0.24. FOLDING FOOTREST*

Stabilising foot support that can also be used as a footboard when moving in and out of the chair.

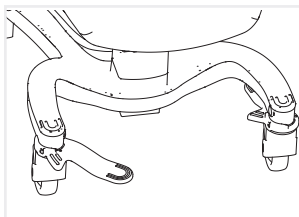


2.0.25. SIDE-MOUNTED SWIVEL FOOTREST*

Mounted on the underside of the front right and/or left chair legs. the footrest can be turned in/out and provide support for a leg as needed. Once the footrest is inside, it is tucked away under the chair frame. The footrest has a rotation stop in both the inner and outer positions.

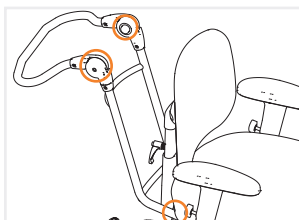

2.0.26. STROLLING BRACKET*

Strolling bracket for transporting a person in the chair. The strolling bracket is mounted to fittings under the seat using two finger screws. This is also where the depth is adjusted. The height is adjusted by pressing the two buttons on the side of the bracket while the handle is positioned at the desired height.

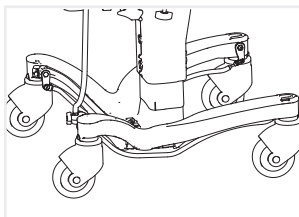

2.0.27. 4-WHEEL BRAKE*

Is activated via the hand/foot brake on all chairs except on chairs with electric brake, where it is activated via hand box or switches. 4-wheel brake brakes on all wheels and provides a stable seat on uneven/skewed surfaces.

Note: On the electric version it is advisable to activate or deactivate the brake completely when using this feature.

2.0.28. ACCESSORIES

This chair is available with a wide range of accessories.
See a complete list of accessories at: www.vela.eu

Note: Charger is included as an accessory when purchasing a VELA chair with electrical functions.

* *Depending on model*

